


## Major II

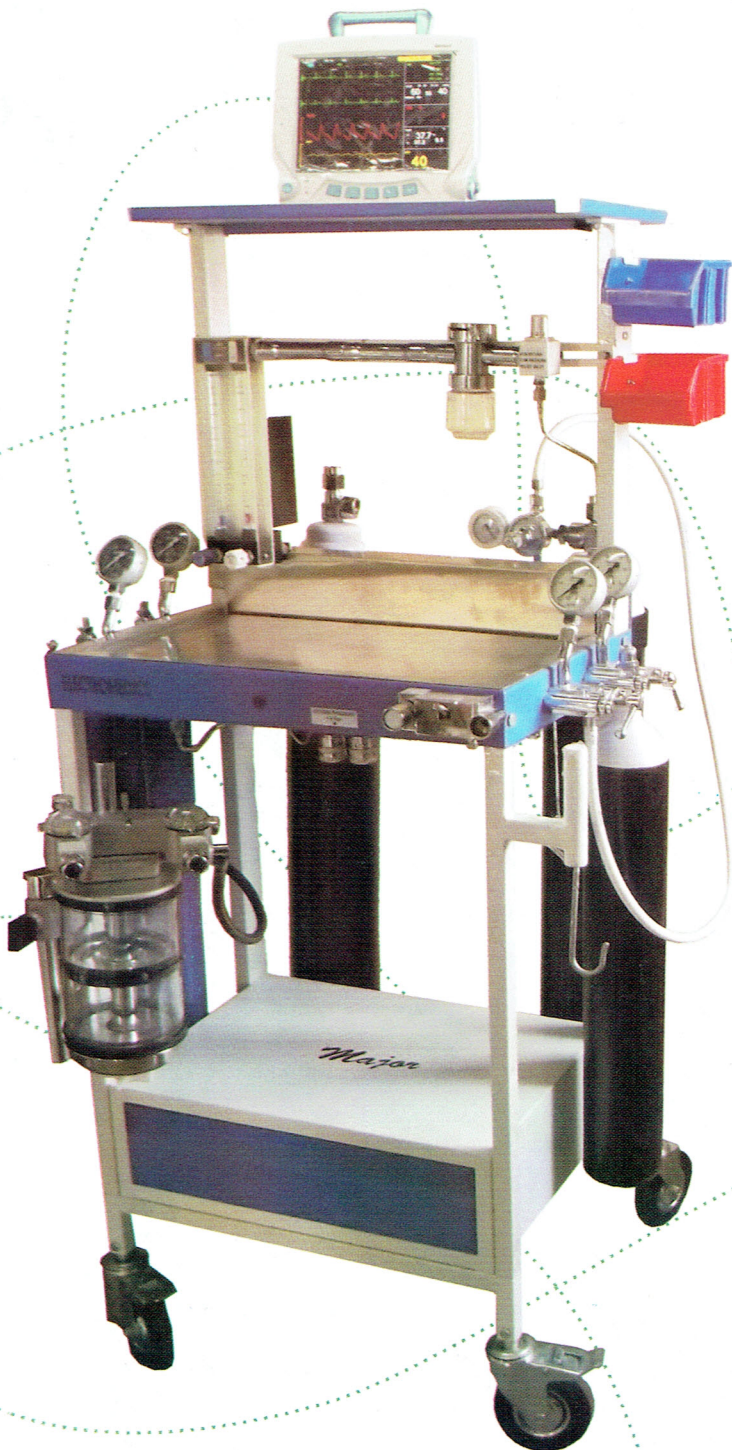

### ◀ Extra Features

- **Construction** : Mild Steel Body, Powder Coated for high corrosion resistance with Anaesthetic appearance.
- **Change Over Unit** : for Change Over open circuit to closed circuit (circle) Emergency Oxygen outlet at table level.
- **Non Return Cum Pressure Relief Valve** : it minimises the risk of backflow of gases and blows off when pressure exceeds 200 cms water column.
- **Pressure gauge** : Superior quality 63mm large diameter.
- **Cylinder Yokes Regulators** : Higher efficiency forged regulators & Yokes with S.S. fitting two for Oxygen and two for Nitrous oxide (six Cylinder carrying Capacity)
- **Table Top Tray** : Stainless Steel Size : 600mm x 355 mm for Anaesthetic Tools & Ampule Bins for injections & medicines.
- **Mobility** : Five 5" diameter anti-static castor wheels with front brakes.
- **Self Sealing Valve** : for Humidifier Bottle or Driving Ventilator. ( Double)
- **Dimension** : H xW x D = 1471mm x 603mm x 585mm.
- **Provision for ( Optional )** : Tek Vaporizer, Hypoxia guard, Circle absorber
- **Provision for** : 2 Extra Vaporizers

Machines are also available in Stainless Steel Trolley

● Authorised Dealer ●

MARKETED BY

**ELECTROMEDICS**

J-4/7, Vighnesh wing, Aditya-Nakoda Enclave, Before Axis Bank, Sinhagad Road, Pune 411 030.

Tel. / Fax : +91-20-24251401 Mobile +91 9850007220 / 90110166 93 / 94 / 95

E-mail : [electromedics@hotmail.com](mailto:electromedics@hotmail.com) | Website : [www.electromedics.in](http://www.electromedics.in)  
Kolhapur - 9011016692